



Follow these steps to download, install and configure MaxTRAQ (2D or Lite+) to record video from a WebCam.

1. Download the MaxTRAQ Installed from Innovision Systems website  
[Download Link](#)
2. Run the Installer to start installing MaxTRAQ
3. Install any software/drivers that came with your WebCam
4. Start MaxTRAQ. If you haven't registered or are using a Dongle, please register your copy of MaxTRAQ to receive your license. Note; You do not need to register if you have a Dongle. You can register from within MaxTRAQ or from [here](#).
5. After MaxTRAQ has started, select File->New Recording.
6. Make sure that Camera Type is set to DirectShow
7. In the Camera Dropdown, select your camera
8. Select the Resolution and the Frame Rate
9. If you're using a Logitech C92X series WebCam, select the following **Advanced Settings** using Camera Settings

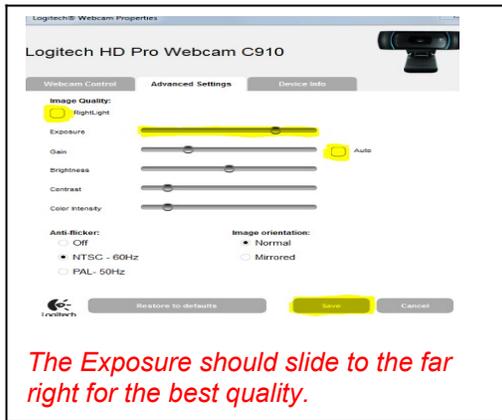

You're now ready to start recording. You can select the number of frames or seconds to record or simply use Record and Stop.

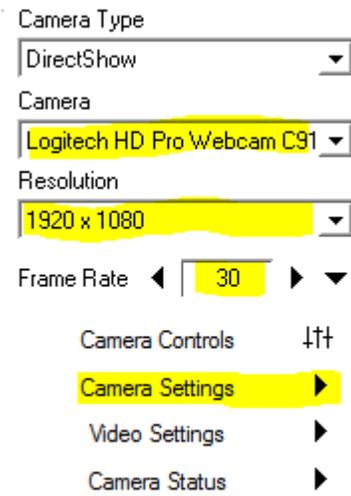